



ONE MORSE IMPLANT

One prosthetic
platform across
all dimensions



CONTENTS

03

ONE MORSE
System overview

04/05

IMPLANT SYSTEM
Connection & product features

06/07

SURGICAL KIT
Instruments & components

08/09

SCIENCE & QUALITY
Surface technology & analysis

10/11

PROSTHETIC WORKFLOWS
Digital & traditional solutions

12

WARRANTY & TITANIUM VIAL
Coverage & implant protection

13


RE-LIFE
Packaging return & recycling

14/15

IMPLANTS & TOOLS
Product range & references

16/17

REGULATORY & LEGAL
Documentation & trademarks

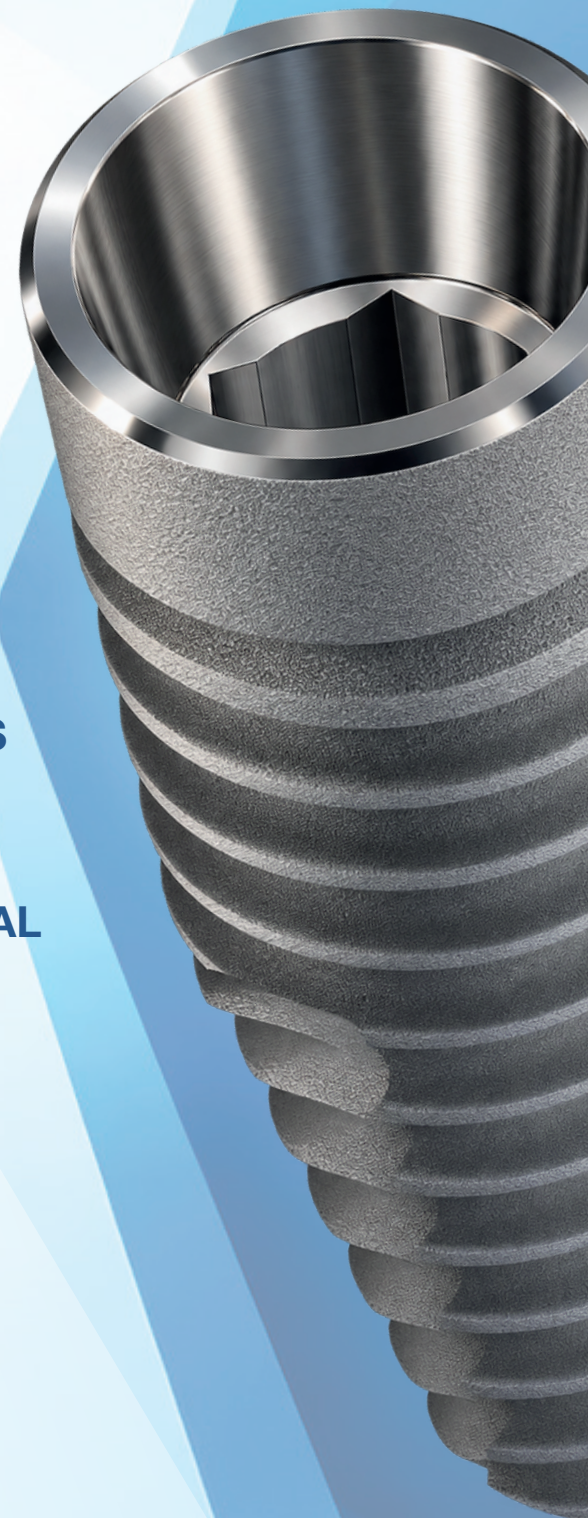

ONE MORSE DESIGNED TO INSPIRE CONFIDENCE

Made from cold-worked Grade 4 titanium, the DESS® ONE MORSE implant combines a tapered, self-cutting design with a spherical apical morphology, promoting high primary stability in lower-density bone and more controlled insertion in higher-density bone, with progressive compaction of the crestal bone.

**DISCOVER
ONE MORSE**
Scan to learn more



One prosthetic
platform



Complete
restorative range



Predictable
results



Reliable
aesthetics

ONE MORSE

CYLINDRICAL PROFILE


Toward the coronal zone, the design becomes more cylindrical with less aggressive threads, allowing smoother insertion, bone compaction, stable seating, and favorable load transfer at the crestal level.

The implant can be placed subcrestal, between 1 and 2 mm, which promotes bone preservation and better long-term integration.

THE CONNECTION

16° Morse taper connection is designed to provide a precise fit between the implant and abutment.

This geometry supports an effective seal and may contribute to long-term mechanical stability.

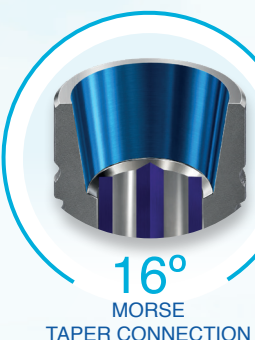



PLATFORM SWITCHING

(crestal load management)

Smaller abutment diameter designed to support platform switching benefits and help optimize crestal load distribution.

Osseointegration Surface Technology

Double treatment: (acid etch & sandblasting) assures optimal bone growth. Complies with SLA standards.

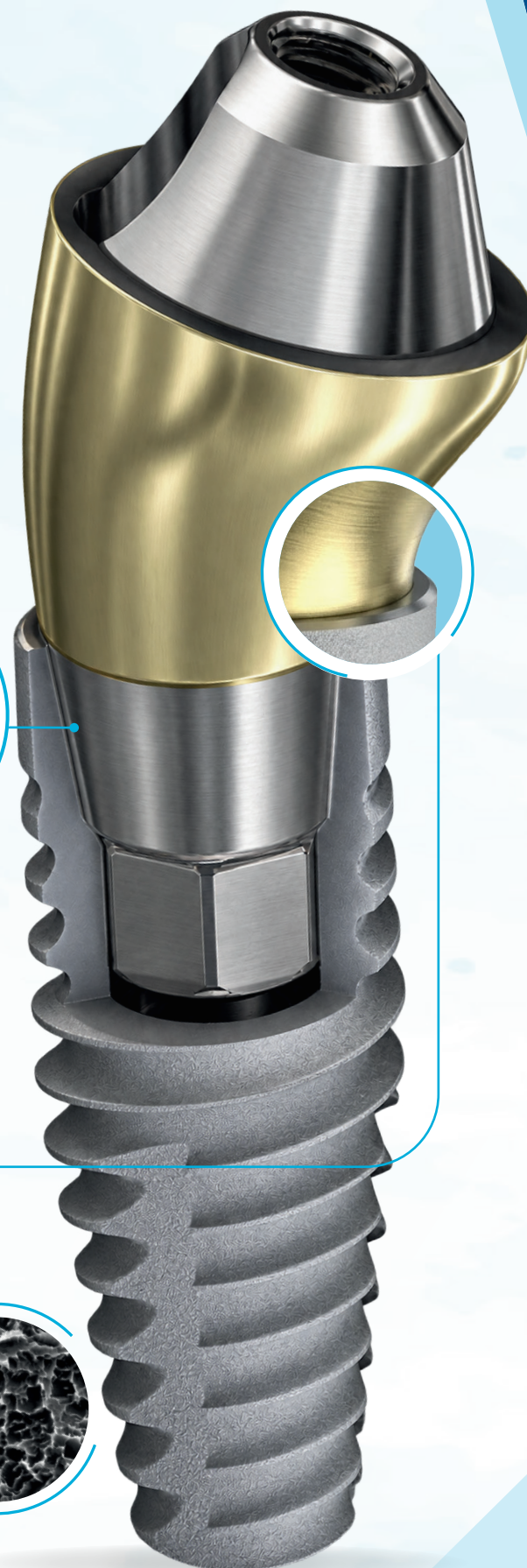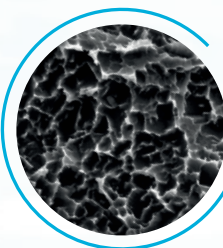


PRODUCT FEATURES

Key connection features supporting fit, stability and long-term restorative performance.

Titanium Grade 4

(mechanical resistance)


- ✓ Cold-worked Titanium Grade 4
- ✓ Provides superior mechanical resistance

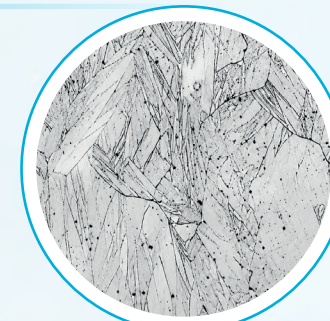

Internal indexing

(guided positioning)


- ✓ Guides positioning during prosthetic placement
- ✓ Helps prevent rotation
- ✓ Supports restorative precision

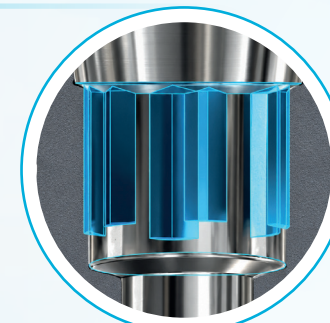

Deep connection geometry

(stable load distribution)

- ✓ Large contact area supports a more even load distribution
- ✓ Designed to enhance long-term connection stability

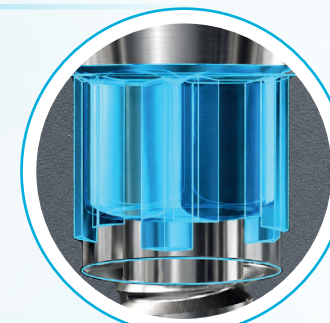

Apical Design

(primary stability)

- ✓ Tapered, self-cutting apical design
- ✓ Supports insertion in soft and very soft bone
- ✓ Helps achieve excellent primary stability



ONE MORSE SURGICAL KIT

STREAMLINED SURGICAL WORKFLOW

The DESS® ONE MORSE surgical kit offers a clear, streamlined set-up for implant preparation, with a modular layout designed to support the selected protocol. Includes the essential instruments for drilling, depth verification and implant placement.

DRILL SIZES

Lengths: 8 / 10 / 11.5 / 13 / 16 / 18 mm
Diameters (color-coded):



Modular & customizable
(adaptable workflow)

- ✓ Arrange the drilling sequence according to protocol and case needs.
- ✓ Supports different clinical indications.



Step-by-step organization
(efficient preparation)

- ✓ Precision, conical and cortical drills with depth gauges and drivers.
- ✓ Clear and intuitive layout for easy instrument identification.



Consistent preparation
(predictable workflow)

- ✓ Streamlined set-up for a safe and efficient surgical procedure.
- ✓ Designed to support predictable clinical outcomes.

SURGICAL KIT CONTENTS

Essential instruments organized for a clear and efficient surgical workflow.

What's included

- ✓ Precision drills
- ✓ Conical drills
- ✓ Cortical drills
- ✓ Depth gauges
- ✓ Implant drivers
- ✓ Screwdriver & handles




IMPLANT AVAILABILITY

	Ø 3.5	Ø 3.75	Ø 4.0	Ø 4.3	Ø 5.0	Ø 6.0	Ø 7.0
Length (mm)							
8.0	✓	✓	✓	✓	✓	✓	✓
10	✓	✓	✓	✓	✓	✓	✓
11.5	✓	✓	✓	✓	✓	✓	✓
13	✓	✓	✓	✓	✓		
16	✓	✓	✓	✓	✓		
18	✓	✓	✓	✓	✓		

THE SCIENCE BEHIND: OST® BY DESS®

Our Osseointegration Surface Technology combines sandblasting and acid etching to create dual-scale roughness.

- ✓ Sandblasting with aluminum oxide particles: for mechanical fixation
- ✓ Double acid etching: microroughness designed to support osseointegration

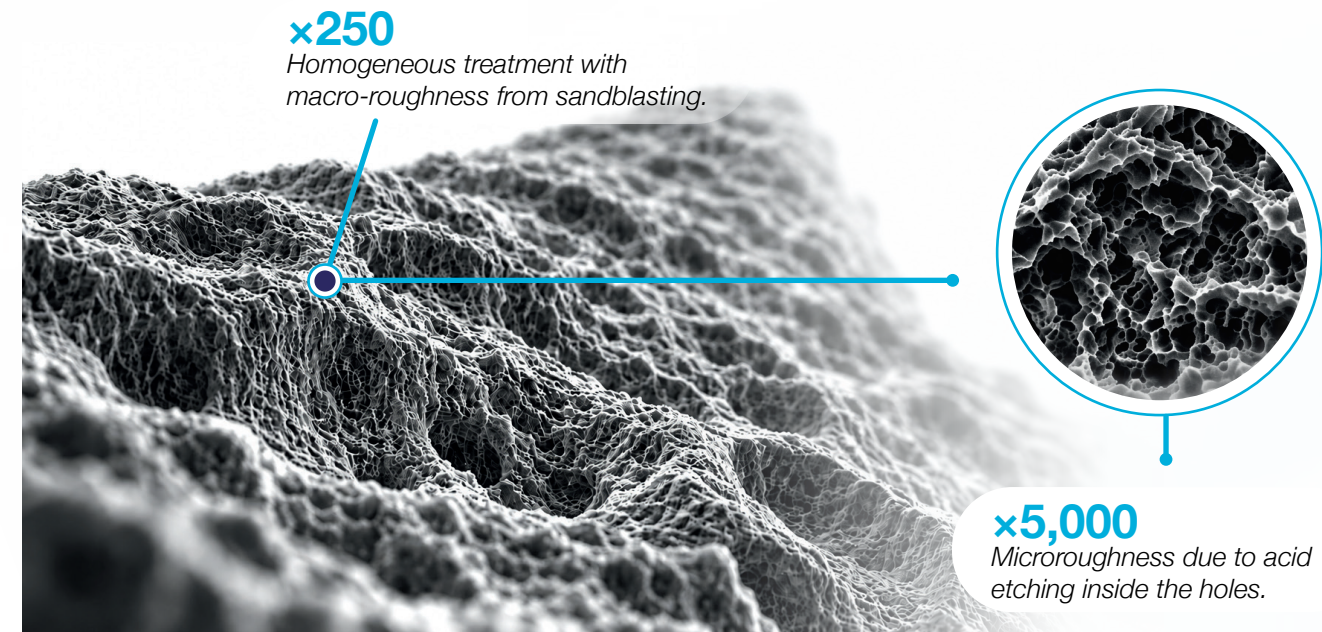



Optimal surface topography for osseointegration.

SURFACE ROUGHNESS BY STEREO SEM

OST® BY DESS®

ACCEPTANCE LIMITS

ISO 25178 parameters

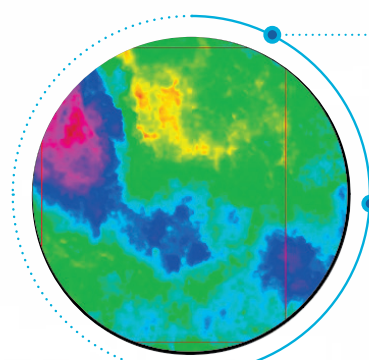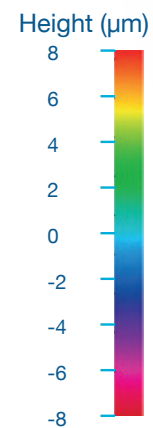

Analyzed area: 60 × 90 µm
Red rectangle: selected area

- Sa**
≥ 1.10 µm
Average height of the selected area
- Sdr**
≥ 30%
Developed interfacial area ratio
- Sz**
> 8 µm
Maximum height of the selected area

RESULTS


The height image obtained by Stereo SEM analysis shows the treated DESS® implant surface. The red rectangle marks the analyzed 60 × 90 µm area.

CONCLUSION

The source report states that treated DESS® implants comply with the defined values for SLA treatment.

ANALYSIS OF SURFACE CHEMISTRY BY XPS

Made from cold-worked Grade 4 titanium, the surface composition analysis shows:

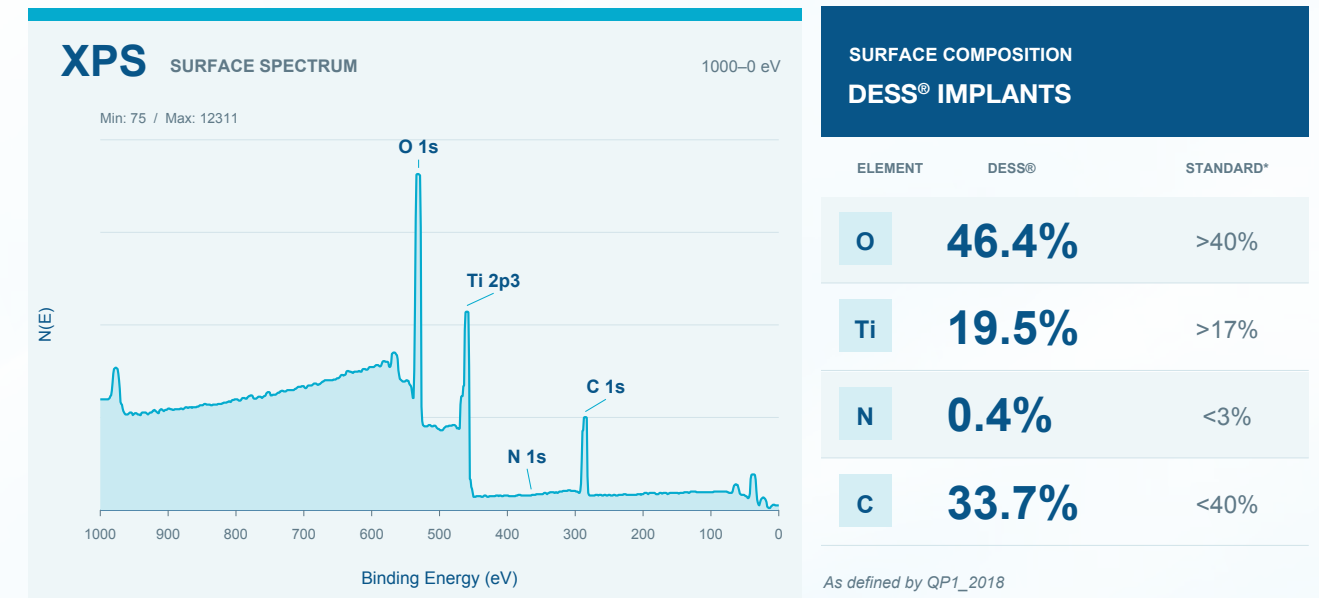

QUALITY & REGULATORY DOCUMENTATION

Our implants are manufactured under ISO 13485:2016 for the design, development and manufacturing of dental implants.

01

ISO 13485:2016

Quality management system



02

CE

Conformity documentation



03

MDSAP


Quality system documentation



04

FDA 510(k)

Reference: K212538

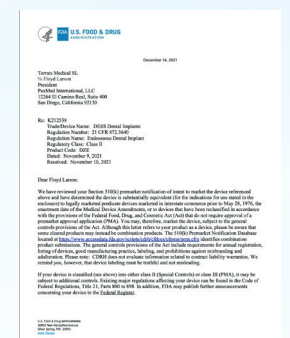

DIGITAL WORKFLOW

01 SINGLE & BRIDGE DIGITAL SOLUTIONS



02 FULL ARCH DIGITAL SOLUTIONS




03 REMOVABLE SOLUTIONS DIGITAL SOLUTIONS



TRADITIONAL WORKFLOW

01 SINGLE & BRIDGE TRADITIONAL SOLUTIONS

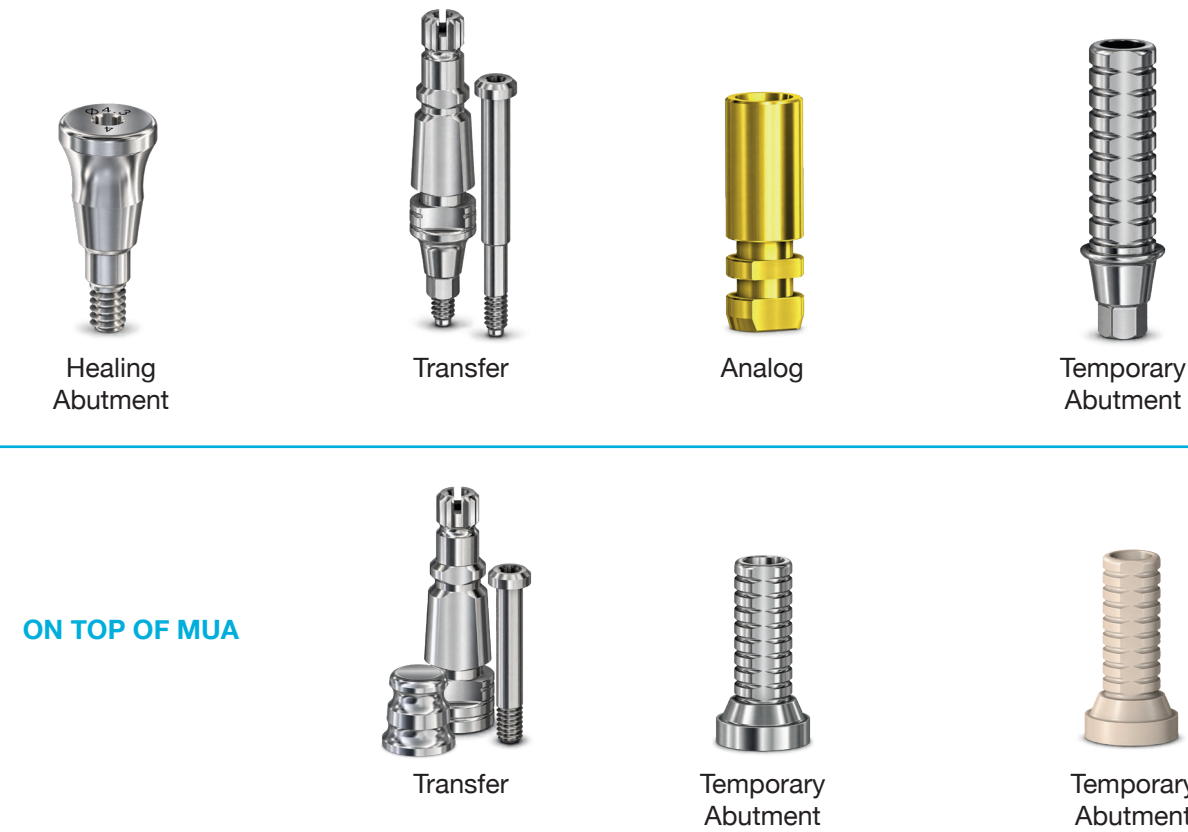


ONE MORSE

One prosthetic platform
across seven implant diameters



DISCOVER
ONE MORSE

Scan to explore
the complete system.

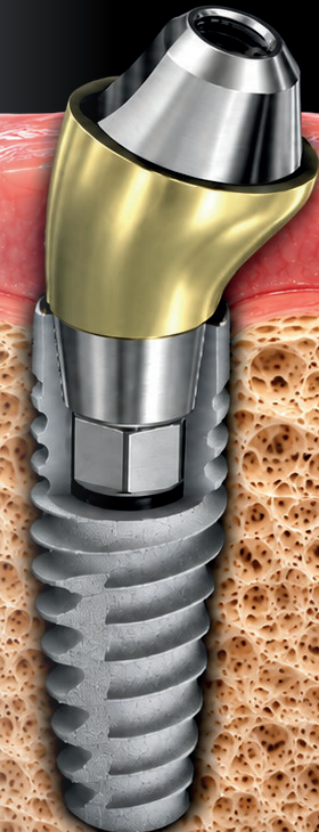

LIFETIME WARRANTY CONFIDENCE IN EVERY COMPONENT

DESS® implants are backed by a lifetime warranty. Coverage also extends to DESS® components and, under specific conditions, to third-party implants used with our components.



WHY CAN WE OFFER A LIFETIME WARRANTY?

- ✓ Rigorous quality controls
- ✓ 100% visual inspection
- ✓ Extensive fatigue testing

Scan to access:



Warranty & Claims



IFUs & Instructions

Important

Third-party implant coverage does not apply once the original implant warranty period has expired. Warranty claims are subject to compliance with the recommendations and instructions provided in the relevant IFUs.

PREMIUM TITANIUM VIAL TITANIUM TO TITANIUM

Designed to protect and preserve the implant in a titanium environment.

Key Benefits


- ✓ Titanium-to-titanium environment
- ✓ Improved implant preservation
- ✓ Reduced risk of contamination from other materials
- ✓ Premium protection system



CLEANING REPORT: TOC

Total Organic Carbon analysis evaluates implant cleanliness.

Method: QIMA-0160
Oxidability / Conductimetry

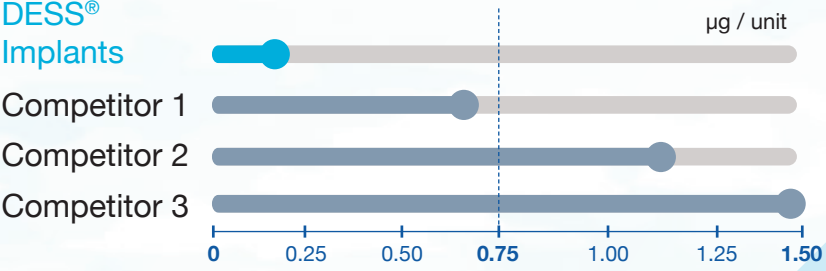

ACCEPTANCE LIMIT ≤ 0.75 µg / unit




GIVE PACKAGING A SECOND LIFE

Return complete implant packaging through Re-Life. DESS® manages collection and recycling.

RETURN THE COMPLETE PACK

Each returned pack should include:

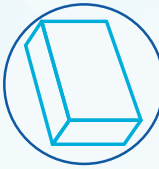

Cardboard box



PET blister



Titanium vial

- 1 RETURN**
Place the empty titanium vial into the blister.
- 2 COLLECT**
Store the complete packs.
- 3 REPACK**
Put the blister back into the cardboard box.
- 4 SEND BACK**
Return them through your usual supplier or distributor.



10 = 1

Return 10 complete packs.
Receive 1 free implant.

RETURN. RECYCLE. REWARD
Together, we can close the loop.



ONE MORSE

IMPLANTS & TOOLS

- Included in the surgical kit Ø = diameter L = length

IMPLANTS

Length	Ø 3.5	Ø 3.75	Ø 4.0	Ø 4.3	Ø 5.0	Ø 6.0	Ø 7.0
8.0 mm	ING0350080	ING0375080	ING0400080	ING0430080	ING0500080	ING0600080	ING0700080
10 mm	ING0350100	ING0375100	ING0400100	ING0430100	ING0500100	ING0600100	ING0700100
11.5 mm	ING0350115	ING0375115	ING0400115	ING0430115	ING0500115	ING0600115	ING0700115
13 mm	ING0350130	ING0375130	ING0400130	ING0430130	ING0500130	—	—
16 mm	ING0350160	ING0375160	ING0400160	ING0430160	ING0500160	—	—
18 mm	ING0350180	ING0375180	ING0400180	ING0430180	ING0500180	—	—

How to read the code: ING 0430 115 = Ø 4.3 × 11.5 mm

DRILLS


L 35 mm	Ø 2.0	Ø 3.5	Ø 3.75	Ø 4.0	Ø 4.3	Ø 5.0	Ø 6.0	Ø 7.0
• Conical	TDRNEO20L	TDRNEO35	TDRNEO37	TDRNEO40	TDRNEO43	TDRNEO50	TDRNEO60	TDRNEO70
• Conical +	TDRNEO20H	TDRNEO35L	TDRNEO37L	TDRNEO40L	TDRNEO43L	TDRNEO50L	—	—
• Cortical	—	TDRNEO35C	TDRNEO37C	TDRNEO40C	TDRNEO43C	TDRNEO50C	—	—



CONICAL DRILL
TDRNEO + Ø




CONICAL DRILL +
TDRNEO + Ø + L

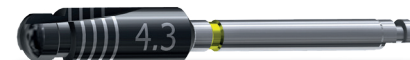


CORTICAL DRILL
TDRNEO + Ø + C

PARALLELIZERS

	Ø 3.5	Ø 3.75	Ø 4.0	Ø 4.3	Ø 5.0	Angled 17°	Angled 30°
•	TPINNEO35	TPINNEO375	TPINNEO40	TPINNEO43	TPINNEO50	TPINNEO20A17	TPINNEO20A30

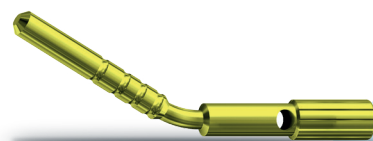

IMPLANT DRIVERS

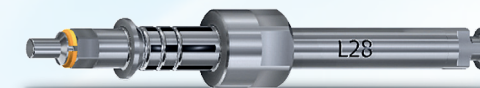

Length

28 mm

37 mm

Ref. 1

• TTRNEOS

• TTRNEOL

Ref. 2

TTRTC28R

TTRTC37R




INSTRUMENTS


- Torque Wrench Tool
DTSTTWC



- Drill Extension
DRILLEXTEND

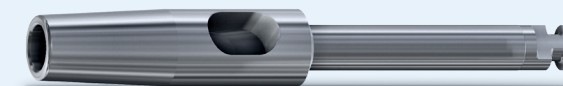

- Special Tools
DTSTCAC
DTSTTRC

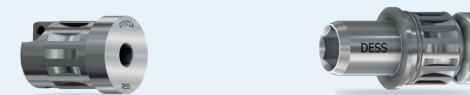


Screw Drivers

- 24 mm DT24NG
- 30 mm DT30NG
- 35 mm DT35NG



COVER SCREW

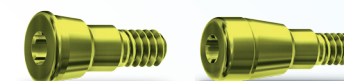

Height

0 mm

20.135-0/0

2 mm

20.135-0/2

HOW TO ORDER

- 1 Choose diameter and length in the implant table (page 14).
- 2 Check the drills and instruments: items marked • are already in your surgical kit.
- 3 Order with the reference code from your DESS® distributor or at www.dessdental.com



All DESS® products are manufactured in accordance with ISO 13485 and are CE marked (Class IIa and IIb medical devices).

All DESS® ANGLEBase® screws and screwdrivers are based on the standard ISO 10664

3i® and 3i Certain®, Miniplant®, MicroMiniplant™ and Osseotite® are registered trademarks of 3i Implant Innovations, USA

ASTRA TECH™, OsseoSpeed™, ANKYLOS® C/X, XiVE®, FRIADENT®, ASTRATECH implant system™ EV, UniAbutment® and Cerec® are registered trademarks of DENTSPLY SIRONA Implants

BioHorizons® is a registered trademark of BioHorizons, INC.

Brånemark System®, NobelReplace® and Replace Select™, NobelActive™, NobelReplace® CC, NobelSpeedy™ and Multi-unit® are registered trademarks of Nobel Biocare AB,

Neodent®, Grand Morse®, Straumann® and synOcta® are registered trademarks of Straumann Holding AG, Switzerland

MIS® is a registered trademark of MIS Implant Technologies Ltd.

Zimmer Screw-Vent® and SwissPlus® are registered trademarks of Zimmer Dental INC, USA

Dyna® is a registered trademark of Dyna Dental Engineering BV

Locator® is a registered trademark of ZEST IP Holdings

Camlog® and Conelog® are registered trademarks of CAMLOG Biotechnologies AG

Megagen® is a registered trademark of MEGAGEN IMPLANT CO. LTD

Torx is a trademark of ACUMENT INTELLECTUAL PROPERTIES, LLC

Bego Semados® is a registered trademark of BEGO GmbH & Co. KG

Anthogyr® is a registered trademark of Anthogyr SAS

Medentis ICX® is a registered trademark of medentis medical GmbH

Osstem®, Hiossen® are registered trademarks of OSSTEM IMPLANT CO. LTD

Dentium® is a trademark of Dentium Co. Ltd.

GLOBAL D® is a trademark of Global D SASU.

Biotech Kontakt® is a trademark of Biotech Dental SAS.

BTI® is a trademark of BTI Biotechnology Institute.

PHIBO® is a trademark of PHIBO Dental Solutions S.L.

DESS products can be used only by licensed dentists, physicians and labs.
Please always refer to the instructions for use available on the website
www.dessdental.com
or supplied with the product.

Some products may not be available in your country. Some of the products are cleared for US and Canada.
Please contact your local distributor for more details.



ONEMORSE-EN/2026

The background is an abstract composition of various shades of blue, ranging from light sky blue to deep navy blue. These colors are used in large, overlapping, angular shapes that create a sense of depth and movement. In the center-left area, there is a cluster of four overlapping, irregular shapes. The topmost shape is a very light blue, followed by a slightly darker light blue, then a medium blue, and finally a solid, vibrant blue at the bottom. These shapes have jagged, torn-paper-like edges, giving them a layered, three-dimensional appearance.

www.dessdental.com